



W H A L T I M A



FOOT AND LEG MASSAGER

MODEL: BORA-BORA
USER MANUAL



Thank you for purchasing this product.

Please read this USER MANUAL carefully before usage.

Read the description of precautions necessary for the proper use of the device with particular care and understanding. Keep this USER MANUAL in a safe place for future reference.







CONTENTS

PRECAUTIONS.....	1
THE NAME AND PURPOSE OF EACH PART.....	3
External structure.....	3
Control panel.....	4
Control panel features.....	5
USER MANUAL.....	5
Pre-use preparation.....	6
After using.....	6
CLEANING AND MAINTENANCE.....	7
Maintenance.....	7
PROBLEM SOLVING.....	8
TECHNICAL SPECIFICATION.....	8

PRECAUTIONS




PRECAUTIONS

- ▶ Before use, carefully read the "WARNING" section to ensure proper operation.
- ▶ The following "Precautions" will help you use the device in the correct way.
- ▶ In order to clearly indicate the degree of danger and damage, sample accidents are classified as "WARNING" and "CAUTION", the instructions marked as follows must be obeyed strictly:

 WARNING	Improper operation of the device may cause injury to its users.
 CAUTION	Improper operation of the device may cause injuries to its users and/or damage to the device.
	 Prohibited. The symbol on the left indicates that independent disassembly is not allowed.
	 Mandatory. The symbol on the left means, that the power cord must necessarily be disconnected after use.

WARNING


The device requires careful supervision when used by children, physically or mentally disabled persons, and when used near designated persons.


 MAKE SURE	<p>Before using the device, it is recommended to consult a special doctor for:</p> <ol style="list-style-type: none"> 1. people who have pacemakers or electronic devices (implants), 2. people undergoing treatment, especially those whose well-being is not good, 3. people with malignant tumors, malignant abscesses, heart disease or high blood pressure, 4. people with skin diseases or skin lesions, 5. pregnant women and during menstruation, 6. people suffering from osteoporosis and/or people who have suffered a spinal fracture, 7. people with high fever.
 PROHIBITED	<p>Children with intellectual disabilities should use the device only under supervision.</p> <p>Do not allow children to use or play with the device by themselves.</p> <p>The device cannot be used by more than one person at the same time.</p> <p>Do not damage, grab, excessively bend, pull hard, twist or bundle the power cord.</p> <p>Do not use the device if the power cord is damaged.</p> <p>Do not bring fingers, hair or other body parts near the nozzle to avoid injury to the user or product damage.</p>
 DO NOT DISMANTLE	No one other than service personnel designated by our company should disassemble, maintain or reconstruct this device. In case of any problems, please contact your local dealer immediately.

PRECAUTIONS





WARNING

 MAKE SURE	Do not use the device continuously for more than 15 minutes.
	Do not place heavy objects on the device.
	The device's default run time is 15 minutes.
	Do not cover the device with a blanket or pillow when using it to prevent the user from being burned and causing a fire by igniting the device.
	Do not step on the device to avoid damage.
	Do not touch the control panel with wet hands.
	Do not leave the device outdoors.
	Do not immerse the entire device in water or other liquids when cleaning.
	If the device is moved from a low-temperature room to a warmer place, it is recommended to start it up after 1 hour to avoid malfunctions caused by condensation.
	Do not use this device in conjunction with other therapeutic devices.
	If the device is used by people who are not sensitive to high temperatures, a temperature-sensitive person must supervise it every time.
	Do not use this device for medical purposes.
	Immediately stop using the device if any abnormal situations (such as the appearance of smoke or a specific odor) occur during use.
	Use only a power supply with a rated voltage.
	Before using the device, inspect it, and if any damage is found, have it repaired immediately by the dealer.

 PROHIBITED	Do not use the device less than half an hour after eating.
	Do not sleep while the device is in use.
	Do not use the device after consuming alcohol or if you have any medical conditions.
	Do not use gasoline, caustic agents or pesticide sprays for cleaning.
	Do not place the device near flammable materials.
	Do not smoke while using the device to avoid the risk of fire.

THE NAME AND PURPOSE OF EACH PART

External structure

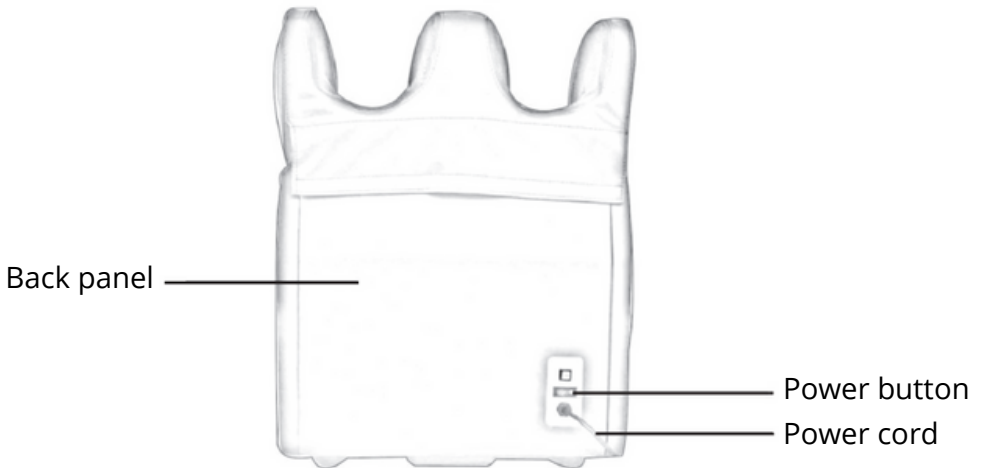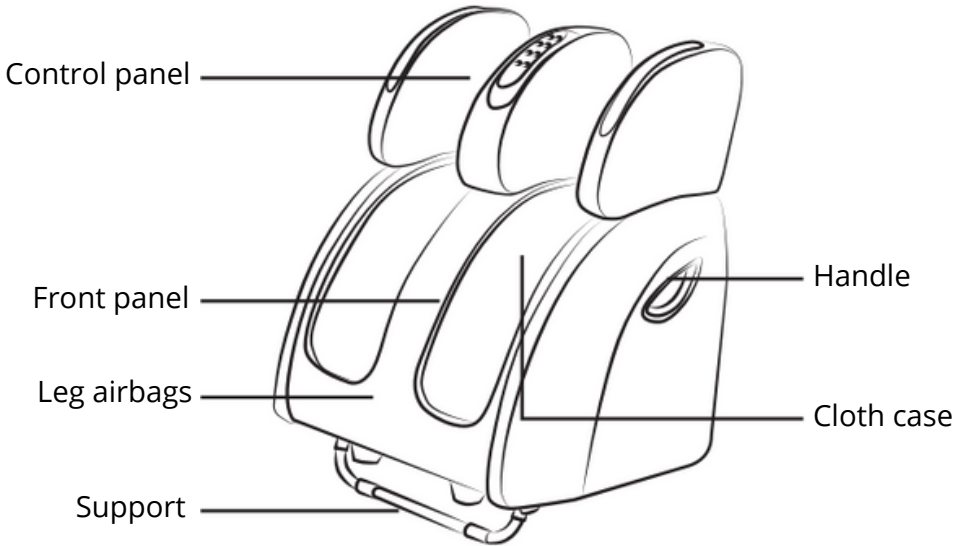

THE NAME AND PURPOSE OF EACH PART

Control panel

Starting at intensity level 2, which can be changed to level 1 or to level 3

Massage program selection for feet

Massage program selection for calves, knees and thighs

Automassage program



Heating function

Vibration function

Rollers designed specifically for sole massage

Power supply switch

INTENSITY



level 1

level 2

level 3



WARNING

- ▶ Do not splash water or other liquids on the control panel.
- ▶ Do not press the button for too long to avoid confusing the program.
- ▶ Do not use the device in a standing position and do not place any heavy objects on it.
- ▶ Do not stand or sit on the device.
- ▶ When transporting the device, do not pull the power cord. Do not pull the unit directly to avoid its damage.

THE NAME AND PURPOSE OF EACH PART

CONTROL PANEL FEATURES

Intensity (Intensity):

Press this button to toggle between medium, high and low intensity settings.

Heating (Heat) button:

If the massager is running, press this button to turn the heating function on/off.

Automatic Foot Massage Mode (Foot):

Press this button to enable/disable the automatic foot massage mode.

Vibration (Vibration):

Press this button to enable/disable the vibration function.

Automatic leg massage mode (Leg):

Press this button to enable/disable the automatic calf, knee and hip massage mode.

Power on/off:

Press this button to turn the power on/off.

Massaging Rollers (Roller):

Press this button to turn the massage rollers on/off.

Auto combination (Auto):

Press this button to enable/disable the auto combination function (no vibration or heating).

Intensity:

Low, Medium, High

Attention:

Medium intensity is the default setting, press the intensity button to switch between 3 settings.

Select the required intensity according to your preference.



- ▶ If a power failure occurs during the massage, turn off the power and unplug the device from the power outlet.
- ▶ Do not bring your hands or fingers close to the unit or foot space while the device is in operation to avoid injuries.
- ▶ Do not let children or pets play with the device to avoid accidents.

Removing the fabric top case:

1. Unpack and remove the fabric case.
2. Clean the fabric case.
3. Put on the top fabric case

(note that the left part of the fabric case is different).

Note that the left part of the fabric case is different from the right.


Removing the lower fabric case:

1. Remove the case with the help of external velcro straps.
2. Unzip the zipper on the case.
3. Remove the fabric case.

Installation of the lower fabric case:


1. Connect the inner velcro straps.
2. Fasten the zipper.
3. Connect the outer velcro straps.

Zipper

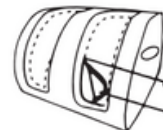

Zipper

Inner velcro strap

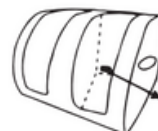


Outer velcro strap

USER MANUAL

PREPARATION BEFORE USE

- ▶ Running the massager:
- ▶ Store the massager in a dry place, away from sunlight.
- ▶ We recommend placing a rug or soft mat underneath to reduce noise during use.
- ▶ Connecting the power supply:
- ▶ Connect the power cord to the device and plug it into an electrical socket.
- ▶ Make sure the power cord is not damaged.
- ▶ Before starting operation, invert the lower support.
- ▶ Turn on the power switch on the edge of the back cover, the power light will turn on, the device will be turned on, and the indicator on the control panel will flash.
- ▶ Press the power button on the control panel to start the device's operation.
- ▶ Select the requested program to start the massage.

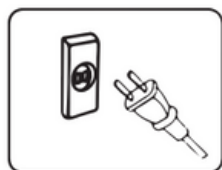

- ▶ Before connecting, check and make sure that the power cord is not damaged.
- ▶ Do not unplug the power cord with a wet hands to avoid electric shock.



Massager scheme 1



Massager scheme 2

AFTER USING

- ▶ Turn the power button on the control panel off.
- ▶ Remove the plug from the electrical socket.

CLEANING AND MAINTENANCE

MAINTENANCE

CLEANING

- ▶ Before cleaning, make sure the cord is not plugged into an electrical socket.
- ▶ Use a soft cloth to clean the surface of the device.



WARNING

- ▶ Never connect or disconnect the plug from the electrical socket with wet hands.
- ▶ Do not use water washers or spray water or other liquids on the inside of the device.
- ▶ Do not use aggressive cleaning agents such as gasoline or thinner to clean the device, otherwise the color of the housing will fade, the device will crack or fade.
- ▶ Do not iron the upholstery.
- ▶ If the device gets wet, do not remove the fabric cover, but let it dry naturally.

MAINTENANCE

- ▶ Make sure the cord is not plugged into an electrical socket.
- ▶ Wipe the dust off. Keep the device clean and away from humid areas.
- ▶ If the device will not be used for a long time, cover it with a dust cover (supplied separately).



WARNING

- ▶ Do not store the device in direct sunlight or in places with high temperatures. This may cause the color of the case to fade.

USAGE TIPS

When using the device on a daily basis, check for the following signs of damage:

- ▶ The smell of burning
- ▶ When the cord is touched, sometimes the power is turned on and sometimes off
- ▶ The power cord is hot
- ▶ Other unusual signs.



WARNING

- ▶ Please stop using the device immediately if any of the signs mentioned above occurs, to avoid damage of the device and possible accidents.
- ▶ Turn the power off and remove the plug from the socket, contact your local dealer or contact us.
- ▶ Do not try to repair or disassemble the device yourself.

FAULT DETECTION

No	Fault symptoms	Possible causes	Possible solutions
1	Failure of startup	The cable is not connected to the electrical socket	Reconnect the plug again
		The power switch has turned off	Turn the switch on
2	The massager does not start and the power indicator is off	The fuse is blown	Insert a toothpick or similar object into the small hole marked "FUSE" above the power switch at the bottom of the rear panel, remove the fuse and replace it with a new one
3	Noise is heard during use	The device has been operating for too long	Turn the power off and wait one hour
4	The external panel is hot	The device has been operating for too long	Turn the massager off and wait half an hour

If any abnormal indications appear after checking the above items, contact your local dealer for assistance from an authorized service engineer.

TECHNICAL SPECIFICATION

Technical parameters:

Device name: Foot and leg massager

Brand: WHALTIMA

Model: Bora-Bora

Rated voltage: 220-240V ~

Rated frequency: 50/60 Hz

Rated power: 40 W

Producer: SHANGHAI RONGTAI HEALTH TECHNOLOGY CORPORATION LIMITED

Location: China, No. 1226 Zhufeng Road, Qingqu, Shanghai

To find the date of manufacture of the device, check the serial number located on the identification label on the device's product case. The serial number consists of 10 characters. The 7th and 8th characters indicate the month, and the 9th and 10th characters indicate the year the device was manufactured.

THE MANUFACTURER HAS THE RIGHT TO CHANGE THE DESIGN, EQUIPMENT AND CHARACTERISTICS OF THE DEVICE IN ORDER TO IMPROVE ITS PRODUCTION.

CHANGES IN TECHNICAL ASPECTS ARE PERMANENT.



W H A L T I M A