W H A L T I M A

**PAMETNI
VAKUUMSKI
MASAŽER**

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Thank you for purchasing this product. Please read this user manual carefully before using to ensure proper installation and operation. Please keep this USER MANUAL properly for future reference.

CONTENTS

01 SAFETY PRECAUTIONS

02 NAME AND FUNCTION OF EACH COMPONENT

Product Structure

03 INSTRUCTIONS FOR USE

Before using

Operation Guide

After using

04 PRODUCT CLEANING AND MAINTENANCE

05 TROUBLE SHOOTING

06 TECHNICAL SPECIFICATIONS

01 SAFETY PRECAUTIONS

- ▶ Before using, please read the SAFETY PRECAUTIONS carefully to ensure correct operation of this appliance.
- ▶ Precautions and warnings listed below will help to ensure proper operation of this appliance. The classes of accidents caused by incorrect use are classified into WARNING and CAUTION.
- ▶ Read all WARNINGS and CAUTIONS below and follow all instructions within this manual.

 WARNING	Improper operation may cause injury to users.
 CAUTION	Improper operation may cause injury to users or damage to product.
 NO Disassemble	 Prohibited

Not suitable for below persons

Close supervision is necessary when this product is used by, on, or near children, invalids, or disabled persons.

 WARNING	<ul style="list-style-type: none">◆ We recommend that those with any of the following conditions should consult your physician before using.◆ Those with pacemaker or electronic devices implanted in the body.◆ Those undergoing treatment, especially those who are not feeling well.◆ Those with malignant tumors, acute diseases, hypertension or heart disease (need to be supervised).◆ Women who are pregnant or menstruating.◆ Those with osteoporosis, distorted joints, dysfunction or recent fracture.◆ Those with serious or injured skin conditions.◆ Those with bad fever.◆ Close supervision is necessary when this product is used by, on, or near children, invalids, or disabled persons.
--	---

Control Button Precautions

 CAUTION	<ul style="list-style-type: none">◆ Do not operate the control buttons with wet hands or pull the adapter wire by force.◆ Do not spill water or other liquid into the control buttons to avoid fault.◆ Do not lay heavy items on the control buttons.
--	---

Best Circumstances

 CAUTION	<ul style="list-style-type: none">◆ Do not use this appliance where temperature is over 40 C .◆ Do not use this appliance in moist or dusty environment.◆ Do not use this product and unplug from outlet during thunder or lightning storms.◆ Do not store near heat or open flame.
--	--

01 SAFETY PRECAUTIONS

Cautions Regarding the Power Connectors and Power Cord

 CAUTION	<ul style="list-style-type: none">◆ Verify that the power voltage stated on the rating label corresponds to your local electricity supply.◆ Do not use a wet hand to unplug the power cord to avoid electric shock.◆ Do not pull on the AC cord, reach to the front plastic end to unplug it.◆ Do not squeeze the power cord by the appliance or other heavy objects.◆ Do not damage, over bend, pull hard or twist the power cord.◆ Please use a matching adapter, please contact the supplier or qualified person for replacement.◆ If the supply cord is damaged, it must be replaced by the manufacturer or its service agent or a similarly qualified person in order to avoid a hazard.
--	---

Best Practices for Safe Use

 CAUTION	<ul style="list-style-type: none">◆ Do not stack heavy items on the product.◆ Do not use other therapeutic device together with this product.◆ Do not put heavy items on the product.◆ Do not use it for prolonged sessions to avoid discomfort. Do not use on the bony area.◆ Do not use the product for therapeutic treatment purpose. Avoid leaving operating product unattended.◆ Do not use the product right after meal, it may cause discomfort.◆ If you feel any abnormal discomfort during the operation of this product, stop use immediately.◆ It is not recommended to use the unit for more than 15 minutes.◆ Check the product spots, if any damages are present, do not use it and refer to a qualified service professional for servicing.◆ If power goes off while in charging, please turn off the product to avoid hazards when the power is back on.◆ Do not use the product on wet skin.◆ This product is not waterproof, Do not spill water or other liquid onto the product.◆ Remove batteries from the appliance before disposing of the appliance.◆ Cut the power supply before removing the batteries.◆ Dispose the batteries properly.
--	---

 PROHIBITED	<ul style="list-style-type: none">◆ Never operate this product if it has damaged/exposed part or power leakage.◆ Do not stand on or stack heavy objects on the product.◆ Do not let a child use or play with the product.◆ Do not fall into sleep while the product is operating.◆ Avoid using the product under the influence of alcohol.◆ Do not put fingers or foreign objects in the openings of the appliance to avoid injury to the user or damage to the appliance.◆ If there is any smoke or unusual smell while using, stop using immediately, unplug the power cord..
 NO DISASSEMBLE	<ul style="list-style-type: none">◆ Do not disassemble this product yourself to avoid product failure.◆ Do not attempt to repair or dismantle the appliance without permission. Please contact us or authorized service agent.
 UNPLUG POWER CORD	<ul style="list-style-type: none">◆ Pull out the power cord from the power strips before maintenance.◆ Never plug in or pull out power cord with wet hand to avoid shock or injury.◆ Pull out the power cord from the power strips while fully charged.◆ Unplug from outlet when the power supply is cut.

01 SAFETY PRECAUTIONS

Best Practices for Safe Use


	<ul style="list-style-type: none">◆ Do not use on the spine, heart, mouth, vagina or soles of the feet,◆ Do not compress the internal organs when using, otherwise accidents may occur or cause physical discomfort. Do not use in the same area for a long time, otherwise it may cause discomfort.◆ Please handle the bag properly to avoid accidents and keep it out of the reach of children.◆ Do not operate under blankets or pillows as overheating may occur, resulting in fire, electric shock or personal injury. Due to the negative pressure gravity will appear scraping phenomenon, it should be postponed for a few days, and then use the instrument, the effect is better. For those with fragile and sensitive skin, please test the oil on the inside of the arm or the root of the ear before using it.◆ Please stop using the oil if it appears red, swollen and itchy. Do not use this device when the body is wet.◆ Keep essential oils in a cool, dry place away from light and out of reach of children.◆ Natural stone needle is thick, need to preheat 5-10 minutes, when the time should pay attention to the temperature of the stone needle after heating.◆ Choose your own suitable temperature range, based on the skin can accept the temperature, do not overheat to prevent burns.
---	--

Best Practices for Safe Battery Use

	<ul style="list-style-type: none">◆ The battery used in this product is a special rechargeable battery and should not be used in other products. Do not charge the removed battery.◆ Do not throw it into the fire to cause excessive heat.◆ Do not Pierce with nails, gravity impact or disintegrate.◆ Do not connect the positive and negative battery terminals directly with metal.◆ Do not store them together with necklaces, hairpins, etc◆ Keep the removed battery out of the baby to prevent swallowing. In case of a mistake, please see a doctor immediately.◆ When battery fluid leaks, do not touch the liquid with bare hands and handle it as follows.◆ When the battery liquid enters the eyes, there is a risk of blindness, Do not rub your eyes with your hands directly. Wash them with water immediately and seek medical attention.◆ When the battery liquid is on the body or clothing, it may cause skin inflammation or injury. Rinse with water. If you still feel uncomfortable, consult a doctor immediately.◆ Only when the product is discarded, it can be decomposed. If you do not comply with it, it may lead to fire, electric shock, injury and other accidents. If you need to repair (except for battery replacement), please consult the store.
---	---

02 NAME AND FUNCTION OF EACH COMPONENT

External structure

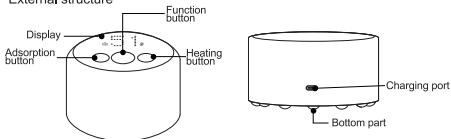Massage diagram

- ▶ The recommended scraping sequence is head, neck, back (thoracic vertebra, lumbar vertebra, sacral vertebra), chest, abdomen, upper limbs, lower limbs (chest, abdomen and other sensitive and fragile parts are not suitable for massage).

- ▶ Each part is generally scraped within 20 times, do not form red spots after the scraping on the skin by force. If the red dots still come out after the scraping, to avoid skin damage, please stop scraping.
- ▶ Scraping strength should be from light to heavy.
- ▶ Usually choose 3 to 5 parts per use is appropriate, when scraping, it should be scraped along one direction.
- ▶ The scraping time should be flexibly adjusted according to each person's age, physical constitution, and the site of the scraping.
- ▶ When using, pay attention to not large area scraping, weak body, only local targeted scraping.
- ▶ Generally, it is necessary to do the second time scraping after the first time scraping 3 to 6 days (the red spots disappear).
- ▶ Normally, there is no strict time and degree restrictions in operation, and it is not necessary to form red spots on the skin by force, 3 to 10 minutes each time, to feel comfortable.

Replace the filter

- ▶ Remove the filter at the bottom of the product to clean or replace.
- ▶ Please pay attention to keep the accessories to prevent loss.
- ▶ The filter can effectively block the essential oil from entering the product during massage, which can effectively extend the life of the product. (After several times of use, it should be cleaned or replaced in time, washed with water, gently wiped with paper towels, and dried in a cool place.)

03 INSTRUCTIONS FOR USE

Before using

- ▶ Before use, please carefully check whether the product is damaged and check whether the accessories are complete.
- ▶ After starting the product, select the product function according to the needs, apply essential oil on the parts of the body that need to be used, and start using.

Battery level display

Low battery	Product display flash
Charging	Negative pressure gear display, inner ring cycle lit up
Fully charged	Negative pressure gear display, inner ring cycle steady on

Operation Guide

On/off

Long press "⊕" for 2 seconds to boot, the screen lights up, showing "01", after starting the negative pressure default 0, the heating default 1. Long press "⊖" for 2 seconds to turn off, and all functions and corresponding displays off.

Adsorption function

In the startup state, press the "⊕" button to adjust the level, 0-1-2-3-4-5-6-0- level cycle is adjustable, the left number of the display screen shows the size of the negative pressure level.

1-2-3, negative pressure mode, negative pressure level increases, "⊕" is on.
4-5-6, breathing mode, negative pressure level increase, "⊕" is on.

Heating function

In the startup state, press the "⊕" button to adjust the heating levels, 0-1-2-3-4-5-6-0- level cycle can be adjusted, the number on the right of the display show the heating levels.

Pressure relief function

In the adsorption state, short press "⊖", the product starts the "pressure relief" function, and in the pressure relief state, short press "⊕" or "⊖", the adsorption gear will return to the pre-pressure relief gear.

After using

- ▶ Turn off the switch when not in use. And place the product safely to prevent damage to the product.
- ▶ When not in use for a long period of time, we recommend that the product packed in the box to keep it away from dust.

CAUTION

- ▶ The product will automatically shut down in 10 minutes, and after 10 minutes, press and hold the switch key again to start up.
- ▶ This product cannot be used when charging, please unplug the charging cable after full charge.
- ▶ Do not operate the controller with wet hands. Otherwise, the controller may fail.
- ▶ Do not splash water or other liquids into the controller to avoid damage to the controller.
- ▶ In order to avoid accidents caused by children accidentally starting the product, the power switch should be turned off when the product is not in use, and then the product should be properly placed.

04 PRODUCT CARE AND STORAGE

Cleaning

- ▶ Be sure to turn off the power and disconnect the power cord from the electrical outlet before cleaning.
- ▶ Use a soft cloth to clean the surface of the product and the controller.
- ▶ To treat stains, spot clean with water-wet or mild neutral detergent-dipped soft cloth, then dry with another cloth.

CAUTION

- ▶ Never use wet hand to plug in or pull out power cord.
- ▶ Never wash this product with water, nor spill water or other liquid into the product.
- ▶ Do not use benzene, thinner or other solvents on the product, it may cause fading or damage to the product.
- ▶ Do not try to dismantle the unit yourself while it's wet, air dry it.

Storage Methods

- ▶ Be sure to cut the power before storage.
- ▶ Place the unit in an environment free of moisture, high temperature and volatile gases.
- ▶ When not in use for a long period of time, we recommend that the product packed in the box to keep it away from dust.
- ▶ When not in use for a long period of time, we recommend that the battery be charged every month to avoid damage.

CAUTION

- ▶ Do not leave this product exposed to direct sunlight for long periods of time. This may cause fading or damage.

Warm Tips

During daily use, please check for the following symptoms:

- ▶ Burnt smell.
- ▶ Power flickers on and off when touching the charging cable.
- ▶ The power cord is hot.
- ▶ Damaged power supply cord.
- ▶ Other abnormal symptoms.

CAUTION

- ▶ Please stop using the product to avoid breakdowns or accidents if any above symptom exists.
- ▶ Switch off the power button and disconnect the plug, contact us for repair.
- ▶ Do not attempt to repair or dismantle the appliance yourself.

05 TROUBLE SHOOTING

No.	Symptoms	Possible Causes	Possible Solutions
1	Noise is heard while in use	The sound is normal	No need to take any measures
2	Appliance does not start	Power Switch is off	Turn on the power switch
		Appliance is overworking and overheat protection cut off the power	Try again after 15mins break
		The battery is dead	Use it after complete charging
		Default timer expires	Restart the appliance

If problem remains, please contact us or authorized service center for maintenance.
Important: Do not attempt to repair this product yourself. Our company shall not be liable for the injury or damage caused by any servicing without authorization.

Disposing the Battery

Batteries can be recycled after use. Please recycle the lithium battery through local programs. Remove batteries from the appliance before disposing of the appliance. Cut the power supply before removing the batteries. Dispose the batteries properly. Battery Removal Steps: Use a professional tool to pry up the top cover piece. Pull out the battery terminals from the main board. Cut the out cover and remove the battery.

06 TECHNICAL SPECIFICATIONS

Name: SMART VACUUM MASSAGER
Model: STONE
Rated voltage: 5V ---
Rated power: 5W
Battery capacity: 2000mAh
Accessories: Type-C cable, user manual

Manufacturer: SHANGHAI AIRONGDA FITNESS TECHNOLOGY CO., LTD.
Address: China, No. 1226 Zhufeng Road, Qingqu, Shanghai.
Made in China.

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Model: STONE

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Producent: SHANGHAI AIRONGDA FITNESS TECHNOLOGY CO., LTD.

Adres: China, No. 1226 Zhufeng Road, Qingqu, Shanghai.

Made in China.

The date of manufacture of the device can be found in the serial number located on the identification sticker on the product body.

THE MANUFACTURER HAS THE RIGHT TO MAKE CHANGES TO THE DESIGN, COMPLETE SET, AS WELL AS IN THE TECHNICAL CHARACTERISTICS OF THE PRODUCT DURING CONTINUOUS IMPROVEMENT OF ITS PRODUCTS.

Do not throw away this manual.

The appearance of the product may not match the image.

This instruction has been thoroughly edited and checked, which, however, does not exclude the presence of errors in the compilation and printing.

We apologize in advance for any possible misunderstanding.

The right of final interpretation belongs to the Company.

This instruction is a literal translation and displays the information stated by the manufacturer.

Improper disposal of waste can harm the environment.

W H A L T I M A